

# REHABILITATION

## Bure Rise & Go Electric Stand Tall Walker



### Bure Rise & Go

The Bure Rise & Go is a cost-effective walker enhanced with a patented power rise function. Thanks to the electric power rise function, Bure Rise & Go becomes a combined stand-up, mobility and walking aid all in one single product.

The Bure Rise & Go enables better work ergonomics for carers, trainers and anyone else working with patients. Thanks to the smart power rise assistance, standing up becomes a safe, confident and natural action that does not wear out the backs of health care personnel.

#### Stable

The walker has anatomical cushions that provide relief for shoulders, arms and the neck and shift the focus to larger muscle groups such as leg and chest muscles instead. This improves stability and makes the walker easier to steer.

#### Easily reachable hand control

The hand controls are located so that both of the carer's and user's hands can easily reach them.

#### Ergonomic handles

Ergonomically designed handles that are steplessly adjustable in every direction for optimum individual adaptation. This is extremely important for stroke patients who can thus reach their "usual" hand position.

#### Folding armrests

Bure Rise & Go's anatomically designed armrests can be folded out to provide the user with an easy means of reaching the handles during the standing manoeuvre.

#### Flexible stepless adjustment

Electrically driven height adjustment makes it easy to pre-set the walker to the desired height. Once the patient is standing in position the height can be steplessly adjusted for comfort and the proper support for an optimal walking configuration.

#### Power Rise

Bure Rise & Go is a unique aid for assisting patients or users during the standing manoeuvre. Personnel burdens, and thus repetitive strain injuries, are minimized.

#### Water proofing

All Bure Rise & Go electrical components have IP 66 international protection ratings and can occasionally be used in damp areas.

#### Stirrup for stability

To provide extra safety Bure Rise & Go is fitted with a stirrup to ensure the walker remains stationary during the standing manoeuvre.

#### Open, reinforced design

The open design provides the same stability but with greater room for manoeuvre. The patient does not feel shut in and it makes assistance with standing up and dressing easier. Bure Rise & Go is also reinforced and extended for extra stability during the standing manoeuvre.

#### Leg support for the standing manoeuvre

Bure Rise & Go is fitted with leg supports that make for a better standing manoeuvre; they are height adjustable and can be moved aside easily to create greater space.

#### Sustainable

The walker is designed to last a long while. For example it has extra reinforcement bands around the arm rest cushions to prevent them from splitting from impacts to the side.

#### Unique power rise harness

The power rise harness is designed to reach as far down across the bottom as possible. This provides optimum support which results in a natural standing action.

#### Lockable castors

Bure Rise & Go has four individually lockable castors as standard.

#### Handbrakes and directionally lockable castors

There are one or two optional handbrakes that can be positioned so that the user can reach them from the walking position. The handbrakes act on the rear castors and the directional locks make it easier for the user to stay on course while walking.

CE tested and approved according to SS-EN 12182:1999

Product Code	Product	Width	Max user weight	Length	Height	Armrests
RHGA-56-316	Bure Rise & Go	77cm	150kg	92cm	97 - 137cm	34 - 67cm



SHA030914

# REHABILITATION

## Bure Rise & Go Accessories



### Bure Rise & Go Accessories

#### Soft Basket

Basket with attachment has enough room for all the important stuff.  
RHGA-56-389



#### Brake Kit Twin

Adds a hand brake so that the user is able to halt movement during walking exercises.  
RHGA-56-334-0




#### Drip Stand with Attachment

Creates great ability for a patient to move freely in a room, on the ward or at home.  
RHGA-56-338



#### Side Support

Prevents arms from sliding from armrest cushions and stabilise walking.  
RHGA-56-337

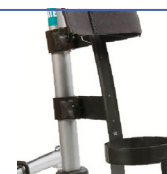

#### Oxygen Bracket

Means that a patient is no longer confined to the bed.  
RHGA-56-336



#### Platform Aluminium - Long/Short

A flexible alternative when moving a patient a short distance. Short or Long.  
RHGA-56-333-1 / RHGA-56-383

Casters with brake  
(various sizes)

#### Walker Gait Belt

Gait Training harnesses are accessories used for increased safety when training using the Bure Rise & Go.  
RHGA-56-388-S/M/L/XL



#### Walker Belt

Washable at 60°, and is simple and flexible to wear, even during visits to the bathroom.  
RHGA-56-384-S/M/L/XL



SHA080415-2